

GTEMP-0225 Rev 1

Replacing the Temple Support Mechanism



Note: You do not have to remove any covers to replace the templesupport mechanism.

1. Description

The temple support mechanism moves the temple supports for the correct positioning of the patient.

Tool Requirements

You must supply the following tools:

- Metric allen keys.
- Metric spanners.
- Flathead screwdriver.



Note: The tool references mentioned in this guide are ISO tool references.

Before Removing or Replacing the Temple Support Mechanism



Before removing or replacing the temple support mechanism, make sure that:

GTEMP-0225 Rev 1


- The unit is disconnected from the main power supply by switching off the circuit breaker in the electrical cabinet.
- The unit is switched off.
- You have the required tools.
- You have the appropriate new temple support mechanism.

2. Removing the Temple Support Mechanism





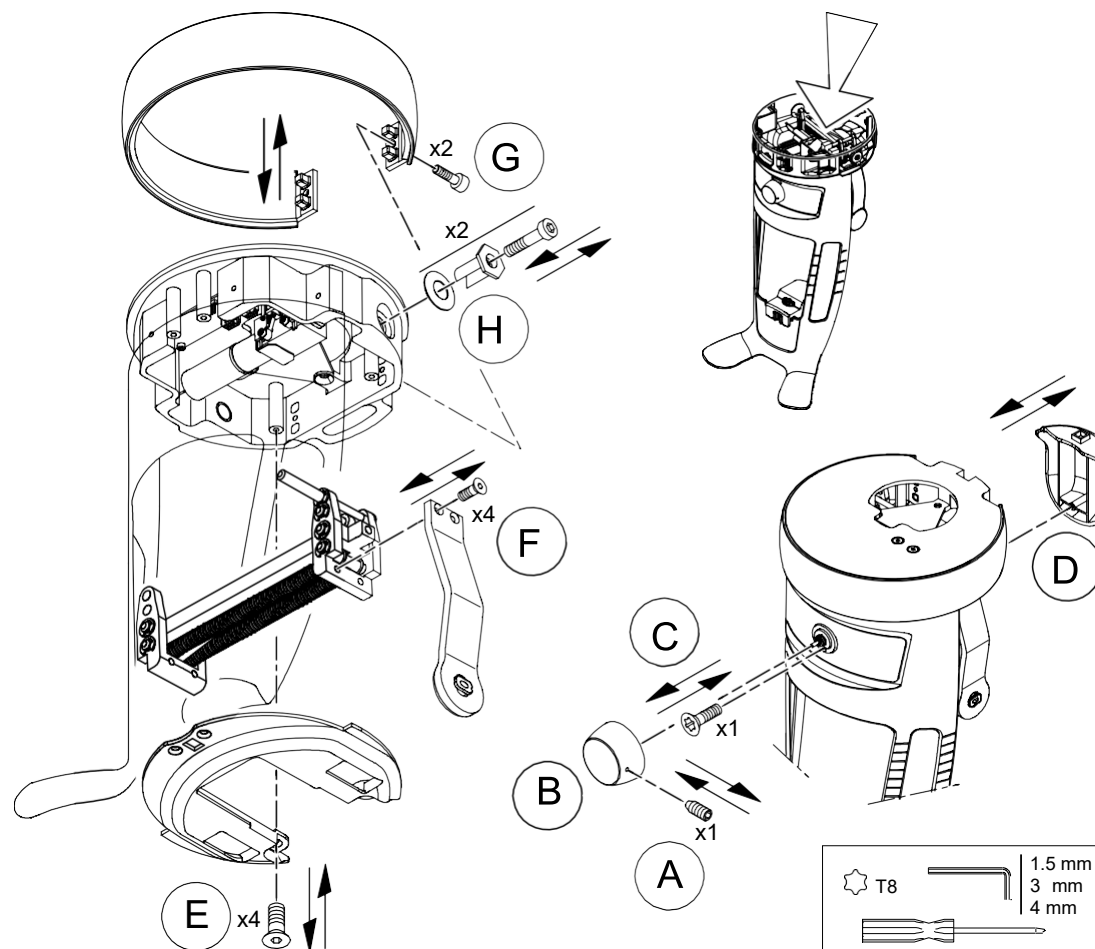

To remove the temple support mechanism, follow these steps:

1. Remove the head and chin rest from the unit.
2. Remove the retaining screw  and the head and chin rest knob. 

GTEMP-0225 Rev 1









GTEMP-0225 Rev 1



3. Remove the screw  and rear cover from the head and chin rest. 
4. Remove the four screws and the head and chin rest cover. 
5. Remove the two screws holding each temple support in place. 
6. Remove the two screws and the head and chin rest collar. 
7. Remove the two screws, sleeves, and washers then remove the temple support mechanism. 

3. Replacing the Temple Support Mechanism

To replace the temple support mechanism, follow these steps:

1. Position the temple support mechanism and insert and tighten the two washers, sleeves, and screws. 
2. Refit the head and chin rest collar and insert and tighten the two screws. 
3. Position the temple supports and insert and tighten the two screws for each support. 
4. Place the head and chin rest cover in position and insert and tighten the four screws. 
5. Place the rear cover on the head and chin rest  and insert and tighten the screw. 

GTEMP-0225 Rev 1

6. Position the head and chin rest knob in position  and insert and tighten the retaining screw. 
7. Replace the head and chin rest on the unit.

4. Testing the Temple Support Mechanism

To test the temple support mechanism, follow these steps using the Service Tools:

1. Run the post installation procedure.
2. Register the intervention.

LAST PAGE OF THE DOCUMENT