



RN-SPCGL47

32

Replacing the Head and Chin Rest



Note: You do **not** have to remove any covers to replace the head and chin rest.

Description

The head and chin rest enables the correct positioning of the patient's head and chin.

Tool Requirements

You must supply the following tools:

- Metric allen keys.
- Metric spanners.
- Flathead screwdriver.




Note: The tool references mentioned in this guide are ISO tool references.

Before Removing or Replacing the Head and Chin Rest

Before removing or replacing the head and chin rest, make sure that:

- The unit is **disconnected from the main power supply** by switching **off** the circuit breaker in the electrical cabinet.
- The unit is switched **off**.
- You have the required tools.
- You have the appropriate new head and chin rest.

Removing the Head and Chin Rest

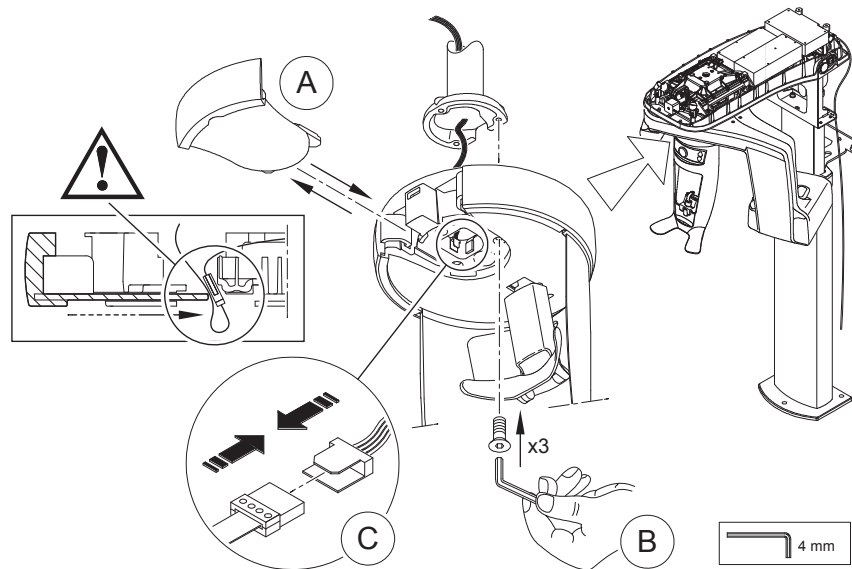

To remove the head and chin rest, follow these steps:

1. Slide off the plastic cover (A).
2. Loosen the three screws (B) but do **not** undo them completely.
3. Disconnect the cable (C).
4. Remove the three screws (B).

Replacing the Head and Chin Rest

To replace the head and chin rest, follow these steps:

1. Slide off the plastic cover (A).
2. Connect the control panel cable (C).
3. Position the head and chin rest on the support bracket.
4. Insert the three screws (B) but do **not** tighten them.
5. Push the control panel cable well into the top of the head and chin rest.

This prevents the control panel cable from getting trapped when you replace the plastic cover.

6. Tighten the three screws (B).
7. Replace the plastic cover (A).

Testing the Head and Chin Rest

To test the head and chin rest, follow these steps using the **Service Tools**:

1. Run the post installation procedure.
2. Register the intervention.
See [“Registering the Intervention”](#) on page 46.