

GTEMP-0225 Rev 1

Replacing the CJ858B Rotation Detector Board



WARNING: Before removing the covers, you **MUST** switch **OFF** the unit **AND** disconnect the unit from the main power supply by switching **OFF** the circuit breaker in the electrical cabinet.

ONLY an approved Carestream technician is qualified to inspect or maintain the unit while it is switched on and while the covers are removed. In this case, **NO** unqualified person must approach the unit.

1. Description

The CJ858B rotation detector board detects the rotation movement.

2. Removing the CJ858B Rotation Detector Board

Tool Requirements

You must supply the following tools:

- Metric allen keys.
- Metric spanners.
- Flathead screwdriver.

GTEMP-0225 Rev 1



Note: The tool references mentioned in this guide are ISO tool references.

Before Removing or Replacing the CJ858B Rotation Detector Board

Before removing or replacing the CJ858B rotation detector board, make sure that:

- The unit is **disconnected from the main power supply** by switching **off** the circuit breaker in the electrical cabinet.
- The unit is switched **off**.
- The head cover is removed.
- You have the required tools.
- You have the appropriate new CJ858B rotation detector board.


Global Replacement Notice

Revision: 01

Title: SPCJ858B Rotation and Y movement detector board Rev 1


Page 3 of 8

GTEMP-0225 Rev 1

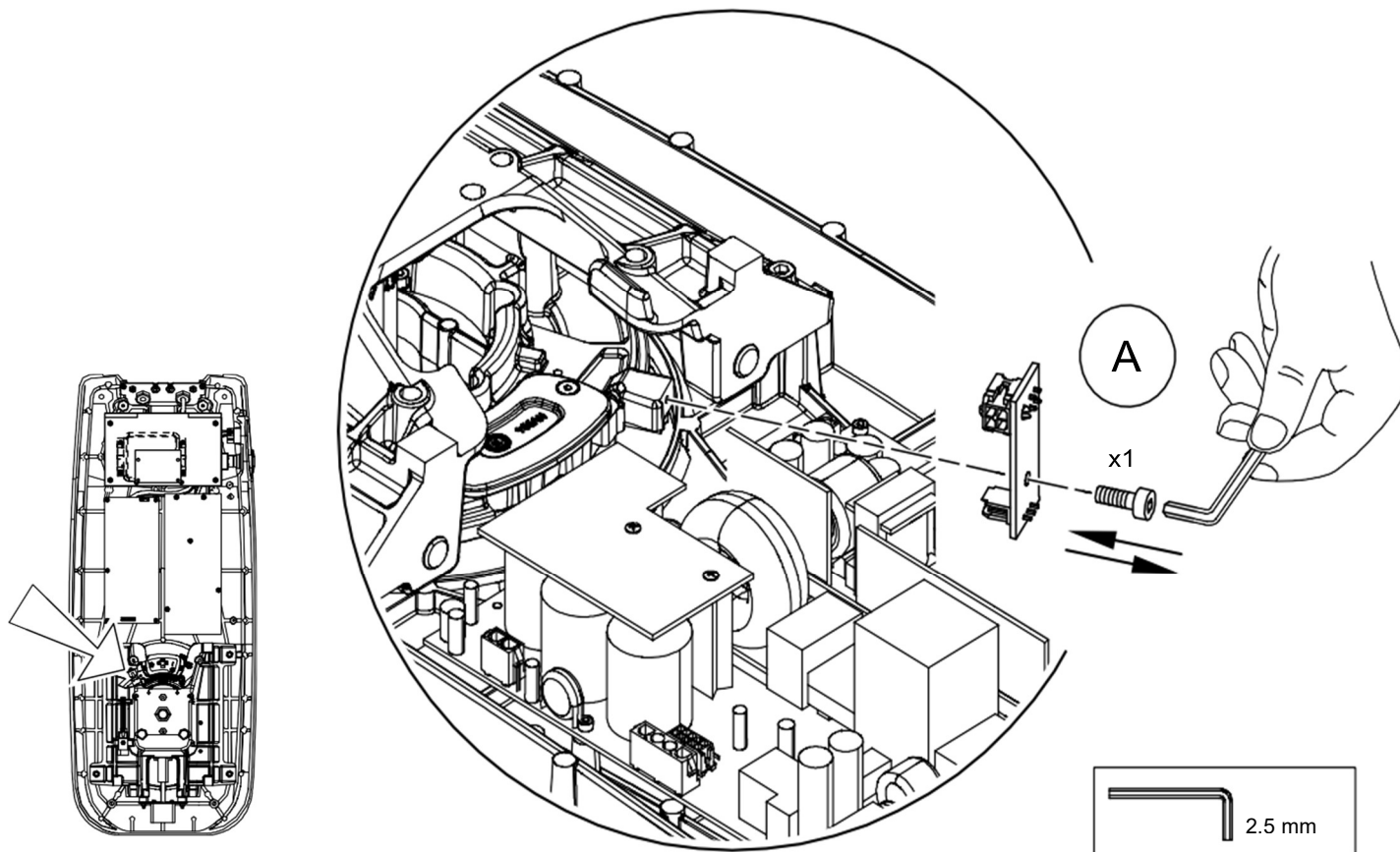


GTEMP-0225 Rev 1

To remove the CJ858B rotation detector board, follow these steps:

1. Remove the screw. 
2. Disconnect the two cables.
3. Remove the CJ858B rotation detector board.

3. Replacing the CJ858B Rotation Detector Board

To replace the CJ858B rotation detector board, follow these steps:

1. Connect the two cables.
2. Position the CJ858B rotation detector board.
3. Insert and tighten the screw. 

4. Testing the CJ858B Rotation Detector Board

To test the CJ858B rotation detector board, follow these steps using the **Service Tools**:

1. Run the post installation procedure.
2. Register the intervention.

LAST PAGE OF THE DOCUMENT

GTEMP-0225 Rev 1

Replacing the CJ858B Y Movement Detector Board



WARNING: Before removing the covers, you **MUST** switch **OFF** the unit **AND** disconnect the unit from the main power supply by switching **OFF** the circuit breaker in the electrical cabinet.

ONLY an approved Carestream technician is qualified to inspect or maintain the unit while it is switched on and while the covers are removed. In this case, **NO** unqualified person must approach the unit.

1. Description

The CJ858B Y movement detector board detects the Y movement.

2. Removing the CJ858B Y Movement Detector Board

Tool Requirements

You must supply the following tools:

- Metric allen keys.
- Metric spanners.
- Flathead screwdriver.

GTEMP-0225 Rev 1



Note: The tool references mentioned in this guide are ISO tool references.

Before Removing or Replacing the CJ858B Y Movement Detector Board

Before removing or replacing the CJ858B Y movement detector board, make sure that:

- The unit is **disconnected from the main power supply** by switching **off** the circuit breaker in the electrical cabinet.
- The unit is switched **off**.
- The head cover is removed.
- You have the required tools.
- You have the appropriate new CJ858B Y movement detector board.


Global Replacement Notice

Revision: 01

Title: SPCJ858B Rotation and Y movement detector board Rev 1


Page 7 of 8

GTEMP-0225 Rev 1

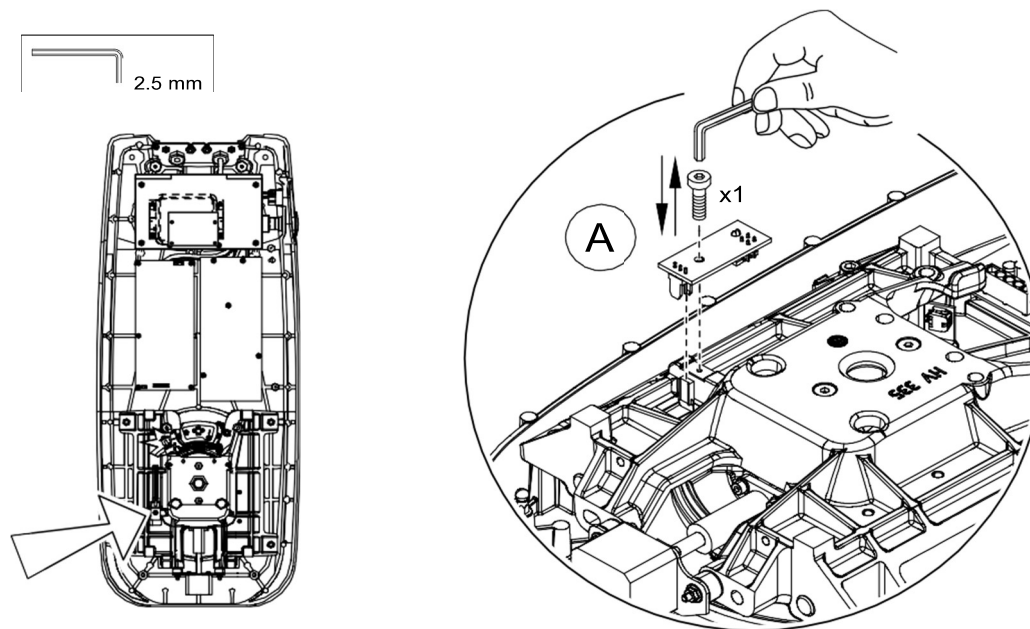


GTEMP-0225 Rev 1

To remove the CJ858B Y movement detector board, follow these steps:

1. Remove the screw. 
2. Disconnect the two cables.
3. Remove the CJ858B Y movement detector board.

3. Replacing the CJ858B Y Movement Detector Board

To replace the CJ858B Y movement detector board, follow these steps:

1. Connect the two cables.
2. Position the CJ858B Y movement detector board.
3. Insert and tighten the screw. 

4. Testing the CJ858B Y Movement Detector Board

To test the CJ858B Y movement detector board, follow these steps using the **Service Tools**:

1. Run the post installation procedure.
2. Register the intervention.

LAST PAGE OF THE DOCUMENT